



The REVIAN RED System includes:

- 1- REVIAN Cap
- 1- Accessory Kit (AC Adaptor, USB Cable, Adjustment Headliners, Travel Bag)
- 1- Charging Support Base
- 1- Downloadable REVIAN app (mobile device not included)

REVIAN specializes in convenient, at-home, light based medical devices that naturally renew your scalp and regrow hair. The REVIAN RED Wireless Cap + Smart App makes hair loss treatment simple, convenient and comfortable in just 10 minutes a day.

Getting Started is as Easy as 1-2-3

1

Charge the Cap with the included magnetic charger for at least an hour to charge and wake it up for the initial setup. A blue light will illuminate around the charging port until fully charged, then softly blink every 4 seconds. The magnetic cable blue light indicates power to the cable and will stay on until unplugged from AC power.



2

Download the REVIAN app with your Android or Apple mobile device. Open the app and click NEW USER to create your account. Accept the requests for bluetooth, notifications, location (Android) and camera during setup.

Follow the instructions in the app to get started. You will register your email address and password with Revian and the system will email you a verification link. Click on the link to accept the verification. Then login to the app using your new credentials.





3

Connect to your Cap from the REVIAN app. To activate, put your finger over the sensor in the forehead area on the inside of the cap (See the image on the right). Open the app and press the yellow Connect button on the treatment page. Once connected, put the cap on your head and press ► to start the treatment.

Important Notes:

1. The Cap must be unplugged from the charging cable before connecting.
2. Turn your mobile device bluetooth on and accept any request to use features in your phone.

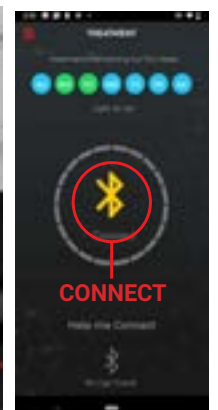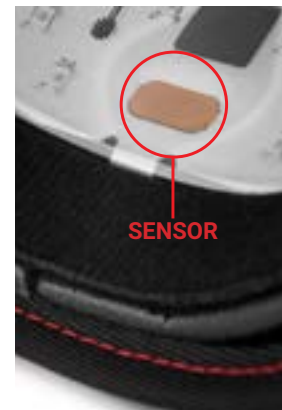


For assistance, please visit us at: REVIAN.com/support or the [Need Help?](#) link in the app slide out tray.

Additional Tips and Recommendations:


Adjusting the Cap may be necessary to ensure a snug fit on your head. A front biosensor awakens the cap and needs to be close against the forehead to operate during treatment. Use the included rear foam pads in the Accessory kit to adjust the fit of the cap.

Pull away the rear foam headliner from the velcro pads and replace with a thicker or thinner foam pad as needed.

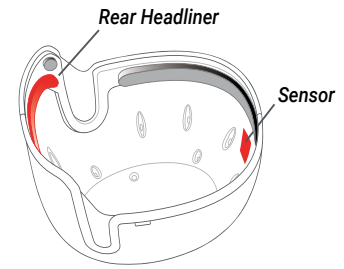

REVIAN App Features

The app has been designed to control your 10-minute daily treatments, guide your progress and help motivate you to stay compliant with your treatments. Here are a few of the features included in the app to help during your hair growth journey.

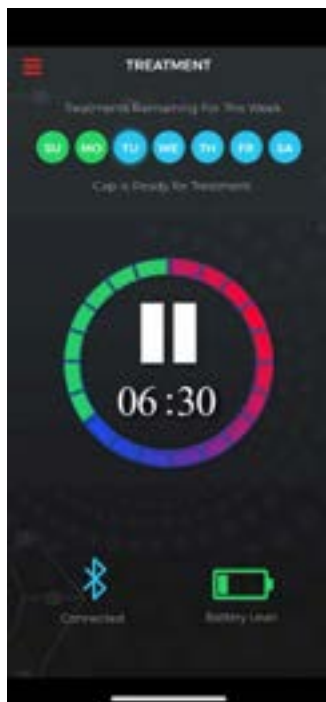

Treatment

- Weekly Treatment Indicators
 - - Missed Treatment
 - - Completed Treatment
 - - Future treatments
- Treatment Play/Pause
- Treatment Timer
- Cap Bluetooth Connection
- Battery Level Indicator

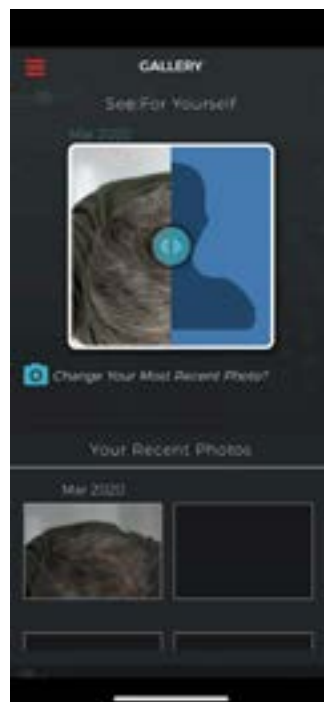


Photo Gallery

- Comparison Slider
- 8 Week Image Gallery
- Past Photos for Viewing
- Retake pictures for 7 days


Seeing is believing so take good quality photos so that you can see the difference over time.



Status

- Motivational Messaging
 - Weekly Tracking
 - Animated Hair Growth Chart
 - Set Daily Reminders
 - Daily Treatment Reminders
- Set treatment reminders at a time that works for you schedule.

Click to View

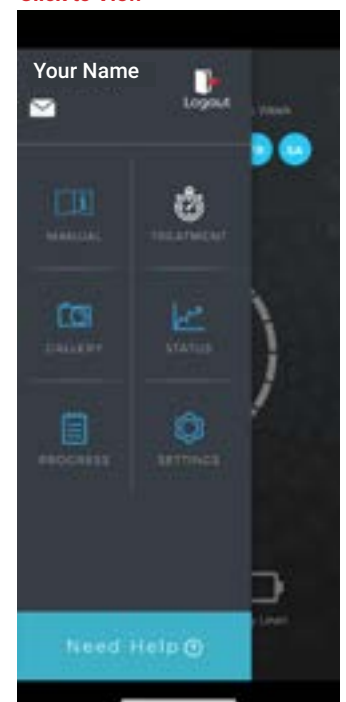

Slide out Tray

- Slide out, Quick Navigation
- Instruction Manual
- Treatment
- Photo Gallery
- Progress Calendar
- Settings
- Logout

For Best Results

- Use your cap for 10 minutes a day. The time of day is up to you, but daily consistency drives good results.
- Take a Day Zero picture (or before picture) of your hair loss area that you most want to improve.
- Set reminders for your daily treatment so you do not miss your treatments.
- Read the full instructions for use in your app or at revian.com/contact/customer-support/